

Linked Page




Certified Register



# Wholesale Supplier

Pediatric (NICU) Equipment (Phototherapy)  
A to Z Hospital Equipment and  
Surgical Products **B to C**

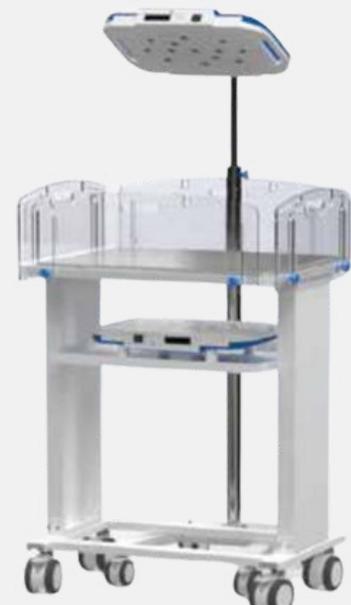

 8275140877

 9822688071

 info@amramedtech.com



Sr No.29/1/2, Adarsh Nagar,  
Uruli Devachi, Hadapsar, Pune  
Branches :- All Maharashtra



amramedtech.com

## STAND MODEL LED PHOTOTHERAPY

### AMAR-SP-0101



18 LED. High power LED provides to cover the bed area and a for effective Phototherapy treatment. No ultraviolet rays ensure to minimize the UV exposures. No infrared rays ensure to minimize the water loss. Adjustable Height - Conveniently allows you to control Radiance levels on the neonate and different allows you to cradle. Convenient Tilting of source Unit, Height Adjustment with Single hand and Castor with elegant maneuvering helps to place the system at different positions, at user convenience. Lamp usage hours greater than 45,000 hours ensure no need for frequent change of lamp.

Weight: 20 Kilograms (kg) (Aprox)

Long Life &gt: 45,000 hrs.

Power LED's emitting blue light in 450-470

HM.

Operation Voltage: 100-250V AC, 50Hz.

Intensity: 14000 LUX

M.S. Fabrication with powder coated paint

## LED DOUBLE SURFACE PHOTOTHERAPY

AMAR-DS-0102



18 LED. High power LED provides to cover the bed area and for effective Phototherapy treatment. No ultraviolet rays ensure to minimize the UV exposures. No infrared rays ensure to minimize the water loss. Adjustable Height - Conveniently allows you to control Radiance levels on the neonate and different allows you to cradle. Convenient Tilting of source Unit, Height Adjustment with Single hand and Castor with elegant maneuvering helps to place the system at different positions, at user convenience. Lamp usage hours greater than 45,000 hours ensure no need for frequent change of lamp.

Weight: 20 Kilograms (kg) (Aprox)

Long Life &gt: 45,000 hrs.

Power LED's emitting blue light in 450-470 NM.

Operation Voltage: 100-250V AC, 50Hz.

Intensity: 14000 LUX

M.S. Fabrication with powder coated paint.

## POWER LED PHOTOTHERAPY



**AMRA-PL-AM-US**



**AMRA-PL-AM-UST**



**AMRA-PL-AM-DSP**

### UPPER SURFACE LED STAND

18 Pc. Hi power LED provides to cover the bed area and for effective Phototherapy treatment No ultraviolet rays ensures minimize the UV exposures. No infrared rays ensures minimize the water loss. Adjustable Height - Conveniently allows you to control Radiance levels on the neonate and different types of cradle. Convenient Tilting of Source Unit, Height Adjustment with Single hand and Castor with elegant maneuvering helps to place the system at different positions, at user convenience. Lamp usages hours greater than 45,000 hours ensures no need of frequent change of lamp.

### DOUBLE SURFACE LED PHOTOTHERAPY

- Optional: Time Totalizer
- Power LED's emitting blue light in 450- 470 NM.
- Long Life > 45,000 hrs.
- Intensity: >25 μw/cm<sup>2</sup>/HM.
- M.S. Fabrication with powder coated paint.
- Operation Voltage: 100-250 V AC, 50Hz

## UPPER SURFACE PHOTOTHERAPY STAND

### AMRA-US-0103

Sturdy mobile stand with height adjustment and easy swivelling of source box  
Unit can be used along with Infant Care Trolley, Radiant Warmer and Incubator.

- Four, 2 Feet imported tubular low pressure mercury vapour fluorescent lamps emitting blue radiation with a peak at 450 nm specially used for hyperbilirubinaemia treatment
- Two day-light tubes for observation
- Clear filter
- Irradiance approx. 12-15  $\mu\text{watts/cm/nm}^2$
- Optional Source Cooling Fan
- Optional Time Totaliser



## UPPER SURFACE PHOTOTHERAPY WITH BABY CRADLE

### AMRA-BC-0104

Sturdy mobile stand with height adjustment and easy swivelling of source box  
Unit can be used along with Infant Care Trolley, Radiant Warmer and Incubator.

Unit can be used along with Infant Care Trolley, Radiant Warmer and Incubator.

- Four, 2 Feet imported tubular low pressure mercury vapour fluorescent lamps emitting blue radiation with a peak at 450 nm specially used for hyperbilirubinaemia treatment
- Two day-light tubes for observation
- Collapsible side supports
- Baby Tray with thick cushion mattress
- Clear filter
- Irradiance approx. 12-15  $\mu\text{watts/cm/nm}^2$
- Optional Source Cooling Fan
- Optional Time Totaliser



## UPPER SURFACE PHOTOTHERAPY WITH TROLLEY

Sturdy mobile stand with height adjustment and easy swivelling of source box  
Unit can be used along with Infant Care Trolley, Radiant Warmer and Incubator.

- Four, 2 Feet imported tubular low pressure mercury vapour fluorescent lamps emitting blue radiation with a peak at 450 nm specially used for hyperbilirubinaemia treatment
- Two day-light tubes for observation
- Clear filter
- Irradiance approx. 12-15  $\mu\text{watts/cm/nm}^2$
- Optional Source Cooling Fan
- Optional Time Totaliser

### INFANT TROLLEY

Sturdy mobile crib mounted on castors can be used along with warmer stand, wall mounted warmer, phototherapy stand & as a separate bed for neonates

Easy access of an Infant from all sides due to clear acrylic collapsible side panels  
Baby tray with transparent slot to incorporate the facility of Under Surface Phototherapy  
Acrylic Baby Tray for infant's head Up / Down facility  
Provision for positioning X-Ray cassette  
Two Lower shelf  
Thick cushion mattress  
Movable on castors

### AMRA-PT-0105



## DOUBLE SURFACE PHOTOTHERAPY UNIT WITH TROLLEY

Sturdy mobile stand with height adjustment and easy swivelling of source box Unit can be used along with Infant Care Trolley, Radiant Warmer and Incubator.

- Four, 2 Feet imported tubular low pressure mercury vapour fluorescent lamps emitting blue radiation with a peak at 450 nm specially used for hyperbilirubinaemia treatment
- Two day-light tubes for observation
- Clear filter
- Irradiance approx. 12-15  $\mu\text{watts/cm/nm}^2$
- Optional Source Cooling Fan
- Optional Time Totalise

### INFANT TROLLEY

Sturdy mobile crib mounted on castors can be used along with warmer stand, wall mounted warmer, phototherapy stand & as a separate bed for neonates

- Easy access of an Infant from all sides due to clear acrylic collapsible side panels
- Baby tray with transparent slot to incorporate the facility of Under Surface Phototherapy
- Acrylic Baby Tray for infant's head Up / Down facility
- Provision for positioning X-Ray cassette
- Two Lower shelf
- Thick cushion mattress
- Movable on castors

### AMRA-DT-0106



**UNDER SURFACE PHOTOTHERAPY** Detachable source box can be used along with radiant warmer and infant crib

Four, 2 Feet imported tubular low pressure mercury vapour fluorescent lamps emitting blue radiation with a peak at 450 nm specially used for hyperbilirubinaemia treatment

Two day-light tubes for observation

Clear filter

Irradiance approx. 12-15  $\mu\text{watts/cm/nm}$

Optional Source Cooling Fan

Optional Time Totaliser

## UNDER SURFACE PHOTOTHERAPY WITH TROLLEY

Detachable source box can be used along with radiant warmer and infant crib

- Four, 2 Feet imported tubular low pressure mercury vapour fluorescent lamps emitting blue radiation with a peak at 450 nm specially used for hyperbilirubinaemia treatment
- Two day-light tubes for observation
- Clear filter
- Irradiance approx. 12-15  $\mu\text{watts/cm}^2/\text{nm}$
- Optional Source Cooling Fan
- Optional Time Totaliser

**AMRA-UT-0107**



### INFANT TROLLEY

Sturdy mobile crib mounted on castors can be used along with warmer stand, wall mounted warmer, phototherapy stand & as a separate bed for neonates

- Easy access of an Infant from all sides due to clear acrylic collapsible side panels
- Baby tray with transparent slot to incorporate the facility of Under Surface Phototherapy
- Acrylic Baby Tray for infant's head Up / Down facility
- Provision for positioning X-Ray cassette
- Two Lower shelf
- Thick cushion mattress
- Movable on castors

## LED PHOTOTHERAPY

An innovative and new generation method, incorporating a state-of-the-art blue LED light source for the treatment of newborn jaundice. Designed to overcome the drawbacks of conventional phototherapy units. Specially arranged LEDs allows to optimal coverage of light over baby.

LED Phototherapy is the next generation LED based technology for treatment of indirect hyperbilirubinemia in newborns. Blue LEDs emit most of their light in the 450-470 nm spectrum. This range corresponds to the peak absorption wavelength (458 nm) at which bilirubin is broken down.

It has features that are needed for intensive, efficacious phototherapy. Advanced blue spectrum High Intensity Super LEDs produce minimal heat with low energy consumption in a dramatically reduced size. Super blue LED's used have long life and low maintenance cost.

- Increased surface area
- High spectral irradiance
- Long lasting blue narrow-band LED light
- Advanced Blue Spectrum Irradiation
- More Efficient Than Other Modes of Therapy
- Compact, light weight, Aesthetic, No heat
- 17 High intensity Super Blue and 4 White LEDs
- Time totalizer with Dose received, Usage, etc

### SPECIFICATIONS

Wavelength	450 nm – 470 nm
Irradiance	< 25 $\mu\text{watts/cm}^2/\text{nm}$
Height Adjustment	Available
Tilting Facility	Available
Power Consumption	25W
Mobile on Castors	

### FEATURES : TREATMENT TIMER

- 2 LINE LCD DISPLAY
- TOTAL USAGE INDICATION
- DOSE TIMING UPTO 24:00H
- MEMBRANE KEYS
- COUNT UP TIMER
- DOSE RECEIVED



### PHOTOTHERAPY UNITS DELUXE


#### LED Bilipod

AMRA-LB-0108




#### ABS

ABS Plastic Construction

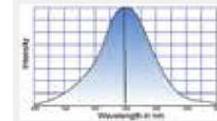

Optimal Wavelength of 460nm

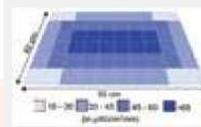

High Irradiance  $>30 \mu\text{W}/\text{cm}^2$  at 45cm distance



No Infrared Radiation



No Ultraviolet Radiation



Low Power Consumption Only 25W



Life of LEDs: 50000 hours



Lamp unit tilt up to 90°



Height Adjustable up to 50 cm



Two Observation Lamp



LCD based timer Patient and Total usage



Fast LED Cooling with Noiseless Fan



Tested as per EN/IEC 60601-1, EN/IEC 60601-1-2 EN/IEC 60601-2-50



European CE Certified


#### Order code:



Bilipod LED Phototherapy Unit From Above



Bilipod LED Phototherapy Unit with Baby Tray

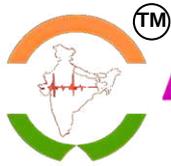

**AMRAMED**<sup>®</sup>  
TECH Pvt. Ltd.

**Wholesale Supplier**


A to Z Hospital Equipment and  
Surgical Products B to C




**Certified Register**

**LED BILIPOD DOUBLE SURFACE DELUXE**

**AMRA-DIB-0109**

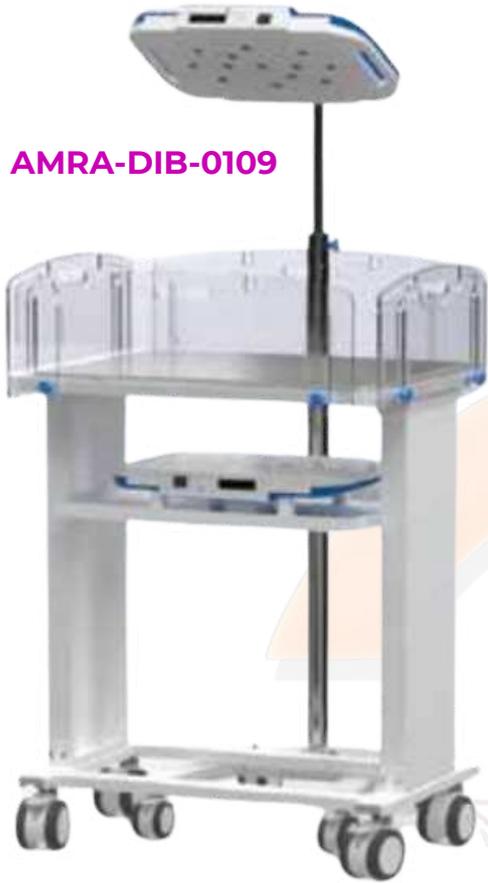

Sturdy Bassinet  
Trolley with  
Dropdown Glasses

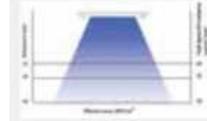

High Irradiance

**Order code:**



Bilipod LED Double  
Surface  
Phototherapy Unit



Available with  
Detachable Trolley

**Order code:**



Bilipod LED Double  
Surface Phototherapy  
Unit with Detachable  
trolley

**Order code:**



Bilipod LED  
Phototherapy Unit  
From Below

**LED BILIPOD 360 DELUXE**


**AMRA-LB-360**



Graphic LCD  
Display



Electrical Elevation  
of Upper couple up  
to 25 cm

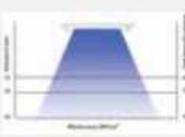

High Irradiance  
>60  $\mu\text{W}/\text{cm}^2$



Slide out tray for  
Easy access of  
the baby



360 degree  
Multi-directional  
Uniform Light  
Distribution



Integrated cooling  
fan system



**AMRAMED**<sup>®</sup>  
TECH Pvt. Ltd.

**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C



**Certified Register**

## LED BILIPOD DELUXE



Winner of



Mother-Baby bonding  
(No barrier during phototherapy)

Portable, Light weight  
(<200g), Space Conserve



Maintenance Free,  
Can be easily cleaned



Power Saving.  
Cost Effective



Easy Installation



Quiet Operation,  
Shock Proof

## LED BILIPOD PORTABLE DOUBLE SURFACE DELUXE

**AMRA-PD-0110**



BILIPAD used for under surface phototherapy

Portable,  
Compact,  
Light Weight



Modern Design,  
Low Space Consumption

Bilipod LED used for overhead phototherapy



Easy Installation



Low Maintenance

Linked Page



Certified Register



**Wholesale Supplier**  
**Pediatric (NICU) Equipment (WARMER)**  
**A to Z Hospital Equipment and**  
**Surgical Products B to C**



 **8275140877**

 **9822688071**


 **info@amramedtech.com**



**Sr No.29/1/2, Adarsh Nagar,  
Uruli Devachi, Hadapsar, Pune**  
**Branches :- All Maharashtra**




**amramedtech.com**

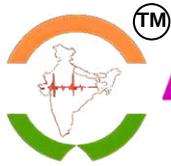


**INFANT RADIANT WARMER**

**AMRA-SW-0201**

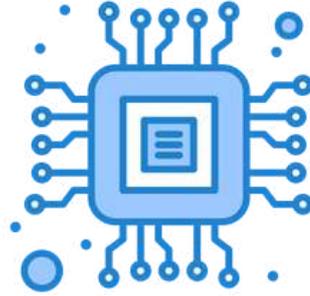


Based on micro processor technology for better accuracy.

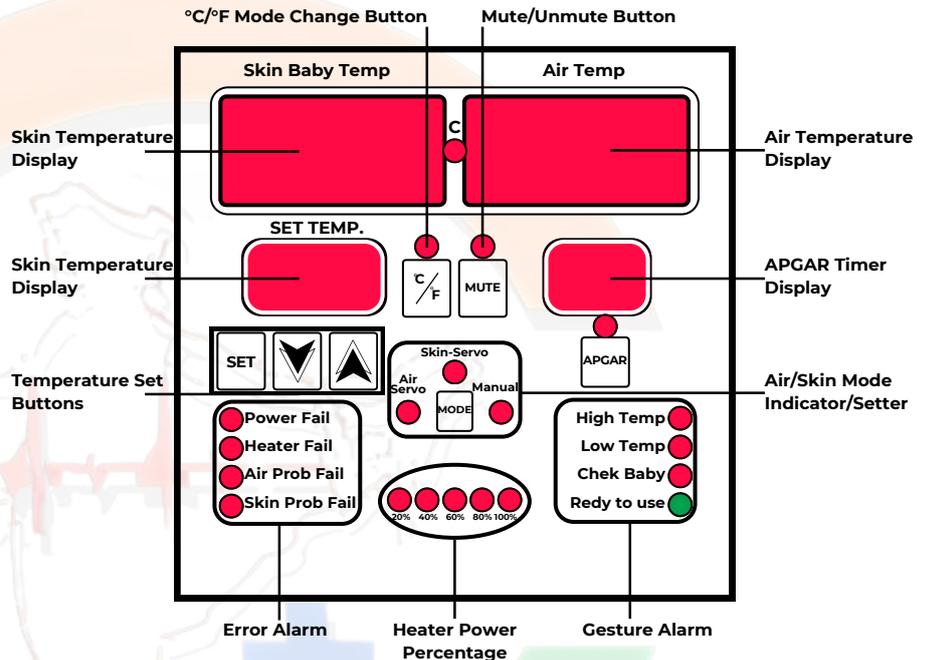

- Three Selection Mode: Serva/Manual/ Prewarm
- Air/Skin Servo Controlled System.
- Acrylic Baby Tray with tilt facilities for infant's head Up/ Down.
- Micro controller technology assure accurate maintenance of set temperature and real temperature.
- Percentage(%) LED Indicator for Heater power
- Heater output with 40% resolution in manual mode
- Alarm indication of Low / High Temprature,Power Failer, Temepreure Prob tailure.
- Mute/ Unmute Facility of alarm.
- Easy access of an infant from all sides due to clear acrylic callapsible side panel.
- Instrument tray, Pole for infusion pump & IV Stand.
- Parabolic refilector to send uniform heat to baby cradle.
- Selt-Hest System.
- User frilendly teather Keypad



**CERAMIC HEATER 220/230V 650W**

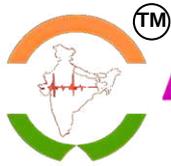

**INTENSIVE CARE SERVO WARMER**

**AMRA-DW-0202**



**AM-WR02(DUAL PROBE) SKIN & AIR PROBE  
SILENT FEATURES**

**EASY TO USE**

As technology advanced, we have kept the Warmer simple & easy to use with touch controls

**FOCUSED HEATING**

The heater with reflector makes it easier to work with the baby by efficiently providing even heat distribution across the cradle and not into the surrounding environment.

**AIR CONTROL MODE**

Air control is an electronic feedback system which functions as a thermostat to maintain a constant temperature of baby cradle by regulating the heat output of Radiant Warmer. It Also acts as a WATCH DOG as Air over temperature cuts off heater even if the Warmer is operating in servo mode

**MANUAL MODE**

The Manual mode allows the Warmer to be kept 'baby ready' without nuisance alarms. When needed the Warmer provides safe, powerful & controlled warming to restore the temperature drop experienced by the baby immediately after delivery.

**TECHNICAL SPECIFICATION**

- Power Source** : 230V AC, +10%,50Hz.
- Light Source** : LED Light to examine babies.
- Heater Capacity** : 600 Watts.
- Power Consumption** : 800 Watts at 100% Heater power.

**SERVO CONTROL MODE**

Servo control is an electronic feedback system which functions as a thermostat a constant temperature at the site of skin probe by regulating the heat output of Radiant Warmer maintaining a constant abdominal Skin temperature between 36.5 & 37.0 C is the simplest way to provide a "Thermoneutral" environment, minimizing the no. of calories needed to maintain normal body temperature & redusing the risks of cold stress or overheating. Servo control is the only acceptable method of heat regulation for the infant cared for under a Radiant Warmer.

**INFANT WARMER (SKIN/MANUAL)**

- Skin temperature displayed on 1 "bright Red LED display".
- Set temperature displayed on separate bright red LED display.
- Two modes of warming i.e. Skin, Air & Manual.
- Heater output in 10 equal steps in percentage.
- Programmable Heater power & Time in Manual mode.
- Individual LED for selection for mode (Skin, Air& Manual).

- APGAR TIMER / ATTEDN BABY ALARM**
- PROGRAMMBLE MUTE TIME SET**
- TEMPERATURE RANGE: SKIN 25C-39C**

- OPTIONAL: X-RAY PROVISION,**
- BATTERY BACKUP, HEIGHT**
- ADJUSTMENT, OVAL SHAPE TROLLEY**



### INTENSIVE CARE WARMER



**AMRA-SW-0203**

**FIXED CRADLE BABY WARMER**

#### AM-WR01 (SINGLE PROBE) SKIN PROBE

#### SILENT FEATURES EASY TO USE

As technology advanced, we have kept the Warmer simple & easy to use with one touch controls

#### FOCUSED HEATING

The heater with reflector makes it easier to work with the baby by efficiently providing even heat distribution across the cradle and not into the surrounding environment

#### SKIN CONTROL MODE

Air control is an electronic feedback system which functions as a thermostat to maintain a constant temperature of baby cradle by regulating the heat output of Radiant Warmer. It Also acts as a WATCH DOG as Air over temperature cuts off heater even if the Warmer is operating in Servo mode.

#### MANUAL MODE

The Manual mode allows the Warmer to be kept 'baby ready' without nuisance alarms. When needed the Warmer provides safe, powerful & controlled warming to restore the temperature drop experienced by the baby immediately after delivery.

AMRA-OH-0204 AMRA-OH-0206  
AMRA-OH-0205 AMRA-OH-0207



**OXYGEN HOOD**

AMRA-LS-0208 AMRA-LS-0209




**LARYNGOSCOPE**

AMRA-AB-0210  
AMRA-AB-0211 AMRA-AB-0212



**SILICON AMBU BAGE**

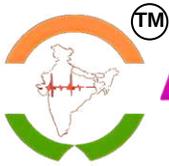

## INTENSIVE CARE SERVO WARMER



**AMRA-SW-0301**

### AM-WR02 (DUAL PROBE) SKIN & AIR PROBE

#### SILENT FEATURES

#### EASY TO USE

As technology advanced, we have kept the Warmer simple & easy to use with one touch controls

#### FOCUSED HEATING

The heater with reflector makes it easier to work with the baby by efficiently providing even heat distribution across the cradle and not into the surrounding environment

#### AIR CONTROL MODE

Air control is an electronic feedback system which functions as a thermostat to maintain a constant temperature of baby cradle by regulating the heat output of Radiant Warmer. It Also acts as a WATCH DOG as Air over temperature cuts off heater even if the Warmer is operating in Servo mode.

#### MANUAL MODE

The Manual mode allows the Warmer to be kept 'baby ready' without nuisance alarms. When needed the Warmer provides safe, powerful & controlled warming to restore the temperature drop experienced by the baby immediately after delivery.

#### TECHNICAL SPECIFICATION

- **Power Source** : 230V AC, ±10%, 50Hz.
- **Light Source** : LED Light to examine babies.
- **Heater Capacity** : 600 Watts.
- **Power Consumption** : 800 Watts at 100% Heater power.

#### SERVO CONTROL MODE

Servo control is an electronic feedback system which functions as a thermostat to maintain a constant temperature at the site of Skin Probe by regulating the heat output of Radiant Warmer maintaining a constant abdominal Skin temperature between 36.5 & 37.0° C is the simplest way to provide a "Thermoneutral" environment, minimizing the no. of calories needed to maintain normal body temperature & reducing the risks of cold stress or overheating. Servo control is the only acceptable method of heat regulation for the infant cared for under a Radiant Warmer.

#### NEANT WARMER (SKIN / MANUAL)

Skin maerature displayed on 1 "bright Red LED display  
Set temperatisplayed on separate bright Red LED display.  
Two modes of Warming Skin. Air & Manual.  
Heater output in 10 equal steps In Pertentage  
Programmable Heater power & Timean Manual mode.  
Individual LED for selection for mode (Skin, Air & Manual)

#### APGAR TIMER / ATTEND BABY ALARM

#### PROGRAMMABLE MUTE TIME SET

#### TEMPERATURE RANGE: SKIN 25° C-39" C

#### OPTIONAL: X-RAY PROVISION, BATTERY

#### BACKUP, HEIGHT ADJUSTMENT, OVAL

#### SHAPE TROLLEY

## INTENSIVE CARE SET-UP

### Photo Therapy Unit (Single Surface)

#### Description

- Fitted with 2 Blue & 4 white fluorescent tubes.(2 feet)
- with baby cradle.
- Mounted on castors.

#### AMRA-ST-0001



### Photo Therapy Unit (Double Surface)

#### Description

- Fitted with 2 Blue & 4 white fluorescent tubes. (2 feet) in upper & lower surface
- with baby cradle.
- Mounted on castors.

#### AMRA-ST-0002



### Radiant Warmer

#### Description

- Skin/Air/Manual(LCD)
- 1"Bright FND, Skin Temp. Display.
- 1"Bright FND, Air Temp. Display.
- LED Bar heater output % Display.
- Set Temp. Display on FND/LCD by mode selection
- All modes Display on LCD.
- Timer Display on LCD.

#### Auto Mode

- All Safety Alarms, Visual display on LCD.


#### Skin Mode

- Skin Temp. High(>0.5°C of setTemp.)
- Skin Temp. Low (<0.5°C of setTemp.)
- Skin Temp. over range (>38°C)
- Skin Temp. under range (<34°C)
- Skin Sensor fail alarm.
- Air Mode / Safety Alarm & High Temp. cut off/ Power failure Alarm.

#### Technical Data:

- Input Voltage 220V/AC+10% 50Hz.
- Light Source 50Watts Bulb to examine Baby

#### AMRA-ST-0004

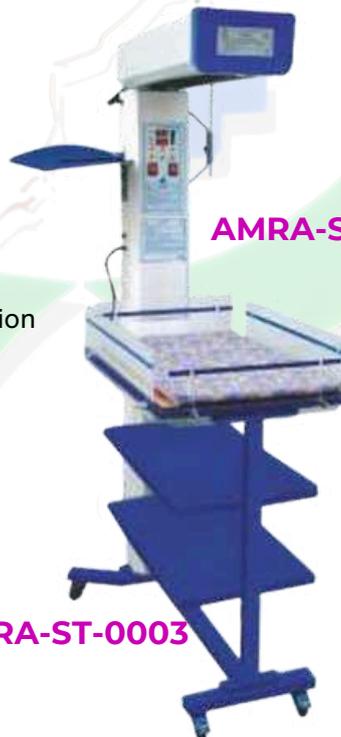

#### AMRA-ST-0005



#### AMRA-ST-0003



**AMRAMED**<sup>®</sup>  
TECH Pvt. Ltd.

**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C



**Certified Register**

## INTENSIVE CARE WARMER

**AMRA-WD-5001**



### **AMRA-WR01(SINGLE PROBE) SKIN PROBE SILENT FEATURES**

#### **EASY TO USE**

As technology advanced, we have kept the Warmer simple & easy to use with touch controls.

#### **FOCUSED HEATING**


The heater with reflector makes it easier to work with the baby by efficiently providing even heat distribution across the cradle and not into the surrounding environment.

#### **SKIN CONTROL MODE**

Air control is an electronic feedback system which functions as a thermostat to maintain a constant temperature of baby cradle by regulating the heat output of Radiant Warmer. It Also acts as a WATCH DOG as Air over temperature cuts off heater even if the Warmer is operating in servo mode.

#### **MANUAL MODE**

The Manual mode allows the Warmer to be kept 'baby ready' without nuisance alarms. When needed the Warmer provides safe, powerful & controlled warming to restore the temperature drop experienced by the baby immediately after delivery

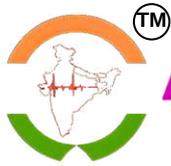

**AMRAMED**<sup>®</sup>  
TECH Pvt. Ltd.

**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C



**Certified Register**

## INFANT WARMERS DELUXE

**AMRA-WD-5002**



10.1" Touch  
Screen



Patient  
Records

01:00:00

Programmable  
Timers



Inbuilt  
Weighing  
Scale



Inbuilt SpO2  
Module  
(Nonin SAT)



Inbuilt  
Resuscitator  
Module



360 degree  
Bed  
Rotation



Smooth Bed  
Tilting



Electronic  
Height  
Adjustment



Polycarbonate  
unbreakable  
Side Support




Modular  
Design



Sliding X-Ray  
Cassette Tray



Molded Fibre  
Heater Unit

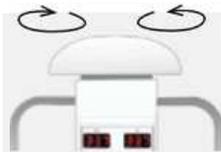

Heater Unit  
rotates 90°  
on both  
sides



LED  
Observation  
Lamp with  
Dimming




Monitor Tray  
Adjusted or  
swivelled



Storage  
Drawers  
provided



Track based  
Rail for Multiple  
Monitors or Pumps

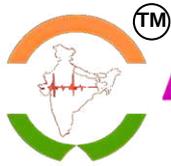

**INFANT WARMERS DELUXE**

**AMRA-WD-5003**



AM IW6000 with  
Neotherm D20  
**AMRA-WD-5004**



AM IW6000 with  
Neotherm D10 Optional




360 degree  
Bed  
Rotation



Smooth Bed  
Tilting



Electronic  
Height  
Adjustment

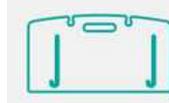

Polycarbonate  
unbreakable  
Side Support




Long Side  
Rail for Multiple  
Monitors or Pumps



Sliding X-Ray  
Cassette Tray



Molded Fibre  
Heater Unit

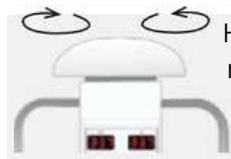


Heater Unit  
rotates 90°  
on both  
sides



LED  
Observation  
Lamp



Monitor Tray  
Adjusted or  
swivelled

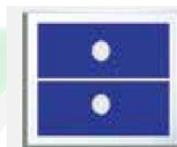

Storage  
Drawers  
provided



European  
CE Certified

Optional



Inbuilt SpO2  
Module  
(Nonin USA)

Optional



Inbuilt  
Resuscitation  
Module

AM IW6000 with Neotherm D20



LED Display  
Skin and Air  
Temp. Control

Order code: AM IW6200

AM IW6000 with Neotherm D10




LED Display  
Skin  
Temp. Control

Order code: AM IW6100

LED  
Phototherapy  
from above



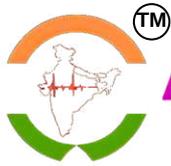

**INFANT WARMERS DELUXE**

**AMRA-PW-6001**



**AMRA-PW-6002**



Optional



Smooth Bed Tilting with X-Ray Tray Facility

**ABS**

ABS Plastic Construction



Heater Unit rotates 90° on both sides



LED Observation Lamp



Polycarbonate unbreakable Side Support



Oval shaped Bassinet Trolley




Lockable Castor wheels



Smooth Bed Tilting



Microcontroller based Temperature Control

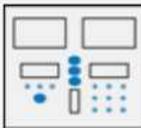

Large Temp Display with Battery Backup



Feather touch User Interface



Battery Backup of Temperature during Power fail




Audio visual Alarms including Sensor fail, Sensor Placement, Heater fail

Optional



LED Phototherapy from above

Optional

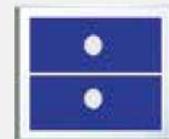

Storage Drawers provided

AM IW7000 with Neotherm D20



LED Display Skin and Air Temp. Control

Order code:AM IW7200

AM IW7000 with Neotherm D10



LED Display Skin Temp. Control

Order code: AM IW7100


**INFANT WARMERS DELUXE**

**AMRA-PW-6003**



**ABS**

ABS Plastic Construction

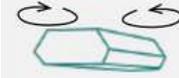

Heater Unit rotates 90° on both sides



LED Observation Lamp



Polycarbonate unbreakable Side Support



Oval shaped Bassinet Trolley



Lockable Castor wheels



Easy Mattress Tilting



Microcontroller based Temperature Control

**AMRA-PW-6004**



Large Temp Display with Battery Backup



Feather touch User Interface



Battery Backup of Temperature during Power fail




Audio visual Alarms including Sensor fail, Sensor Placement, Heater fail

Optional



LED Phototherapy from above

Optional

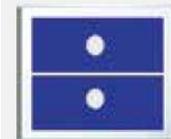

Storage Drawers provided

Optional



Smooth Bed Tilting with X-Ray Tray Facility

AM IW5000 with Neotherm D20



LED Display Skin and Air Temp. Control

Order code: AM IW5200

AM IW5000 with Neotherm D10




LED Display Skin Temp. Control

Order code: AM IW5100

**INFANT WARMERS DELUXE**

**AMRA-HD-7001**

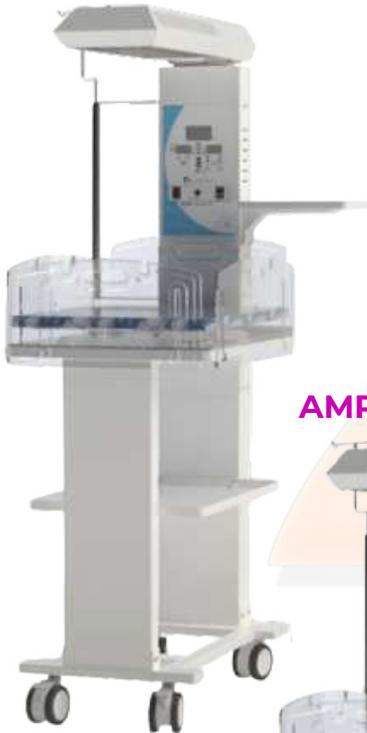

AM IW4000 with Neotherm D10

**AMRA-HD-7002**



AM IW4000 with Neotherm D20

**MS**

MS Plastic Construction



Heater Unit rotates 90° on both sides



LED Observation Lamp




Polycarbonate unbreakable Side Support



Rectangle Bassinet Trolley



Lockable Castor wheels

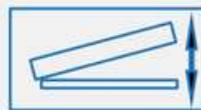

Easy Mattress Tilting



Microcontroller based Temperature Control



Large Temp Display with Battery Backup



Feather touch User Interface



Battery Backup of Temperature during Power fail




Audio visual Alarms including Sensor fail, Sensor Placement, Heater fail

Optional



LED Phototherapy from above

Optional

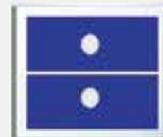

Storage Drawers provided



Tested as per EN/IEC 60601-1, EN/IEC 60601-1-2 EN/IEC 60601-2-21

AM IW4000 with Neotherm D20



LED Display Skin and Air Temp. Control

AM IW4000 with Neotherm D10



LED Display Skin Temp. Control



European CE Certified

Order code: AM IW4200

Order code: AM IW4100



**AMRAMED**<sup>®</sup>  
TECH Pvt. Ltd.

**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C



**Certified Register**

**INFANT WARMERS DELUXE**

**AMRA-HD-7003**



Smooth Bed  
Tilting



Sliding X-Ray  
Cassette Tray

AM IW4000A with Neotherm D20



Fixed  
Bassinet  
Trolley



LED Display  
Skin and Air  
Temp. Control

Order code: IW4200A

AM IW4000A with Neotherm D10



LED Display  
Skin  
Temp. Control

Order code: IW4100A

**AMRA-HD-7004**



Smooth Bed  
Tilting



Sliding X-Ray  
Cassette Tray

AM IW4000AD with Neotherm D20



Detachable  
Bassinet  
Trolley



LED Display  
Skin and Air  
Temp. Control

Order code: IW4200AD

AM IW4000AD with Neotherm D10



LED Display  
Skin  
Temp. Control

Order code: IW4100AD

**INFANT WARMERS DELUXE**

**AMRA-SM-021**  
Stand Warmer



**AMRA-WM-022**  
Wall Mounted Warmer



Order code:

IW4200 S-IW4000 Stand Warmer with  
Neotherm D20 Control

IW4100 S-IW4000 Stand Warmer with  
Neotherm D10 Control

Order code:

IW4200 W-IW4000 Wall  
Mounted Warmer with Neotherm  
D20 Control

IW4100 S-IW4000 Wall Mounted  
Warmer with Neotherm D10  
Control

**INFANT WARMERS - TEMPERATURE CONTROLS**

**Neotherm AM D10**



**34.4** LED based Large  
Skin Temperature  
Display

 Audio visual Alarms  
including Sensor fail.  
Sensor Placement,  
Heater fail

**MODES** Manual  
and Auto


**Neotherm AM D20**



**34.4** LED based Large  
Skin and Air  
Temperature  
Display

 Audio visual Alarms  
including Sensor fail.  
Sensor Placement,  
Heater fail

**MODES** Manual  
Baby Air

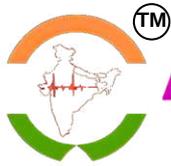

## INFANT WARMERS DELUXE

### Optional Modules

#### AMRA-2 IN 1

#### Inbuilt LED Phototherapy

Six High power LED Lamps for intensive phototherapy treatment are placed near the heater unit.

Can be included with any model of Infant Warmer



#### Inbuilt Resuscitation Module


Option 1:

Easy Resuscitation with Air Oxygen Mixer, Pressure Manometer, and Max Pressure Valve Can be used with T-Piece

Option 2:

Easy Resuscitation with Air Oxygen Blender, Pressure Manometer and Max Pressure Valve Can be used with T-Piece

Slow Suction Unit can be ordered separately

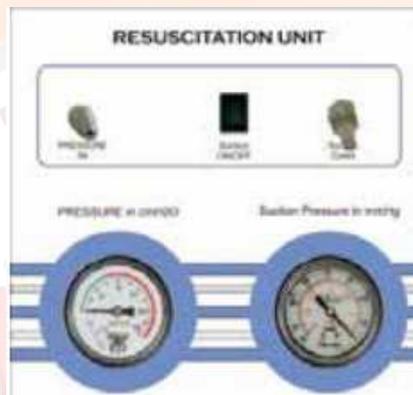

## INFANT WARMERS DELUXE



#### Inbuilt Drawer Box

2 Pull out Drawers made of MS

Can be configured with any model of Warmer

#### Inbuilt PulseOx Module



NONIN PUNISAT  
Technology  
Superior Accuracy in Low Saturation and Dark Skin Patients \*



#### Inbuilt Weighing Scale

Tare Facility given

Accuracy: 10 gm



**PREMATURE BABY INCUBATORS DELUXE**

**AMRA-BI-8001**



Optional

PB11100 with Neotherm 010  
PB11200 with Neotherm D20  
PB11400 with Neotherm iWarm



Smooth Bed  
Tilling

PB11400 Only



5.4" Color Touch  
Screen TFT  
Display

PB13400 Only



Electronic  
Height  
Adjustment



LED based  
Large Skin and  
Air Temperature  
Display

**MODES**

Manual  
Baby Air

**Bassinet Trolley**

**AMRA-BT-8002**



**AMRA-BT-8003**

