


Linked Page



Certified Register

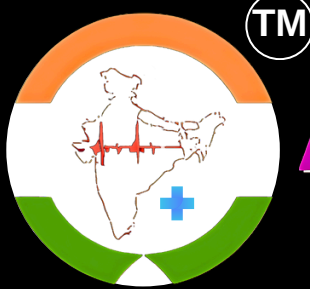

TM

**AMRAMED**<sup>®</sup>  
**TECH Pvt. Ltd.**

**Wholesale Supplier**  
**O.T. TABLE DELUXE**  
**A to Z Hospital Equipment and**  
**Surgical Products B to C**



 **8275140877**

 **9822688071**

 **info@amramedtech.com**



**Sr No.29/1/2, Adarsh Nagar,  
Uruli Devachi, Hadapsar, Pune**  
**Branches :- All Maharashtra**



**amramedtech.com**

## ELECTRIC OT TABLE



### AMRA MATRIX

#### Special Features

- We are providing SMPS based control box which can handle better Voltage Fluctuations, When SMPS is idle there is zero current consumption and saves power. These are as per new European and American standards where the idle power consumption should be zero and where transformer based system are outdated.
- We provide remote which is ergonomically designed to have an easy and better grip with spiral cable and it has LED backlit switch with symbolic position figures making it convenient to use even in the dark.
- We provide zero position facility i.e by pressing the single button OT Table will come to normal position.
- We provide patient reverse orientation mode when head and leg section are interchanged.
- We provide 2 memory function in control panel. • We provide indication of battery back in the second handset.
- The Handset is designed with Back LED Indicator for better visibility.
- Head and foot section is manually operated by the means of Ratchet System.
- We provide Stainless steel Covered Base and Column Covers for easy cleaning and hygiene.
- We provide Anti Bacterial and Anti Fungal Mattress.
- The motion of table is easy with double swivel castors 100mm diameter, in which 2 castors having brakes

#### STEEL QUALITY 304 MEDICAL GRADE STAINLESS STEEL

#### Technical Data

|                                   |                      |
|-----------------------------------|----------------------|
| Length of Table                   | 1900 mm              |
| Width of Table Top                | 520 mm               |
| Minimum Height (without Mattress) | Min. 750 mm          |
| Maximum Height (without Mattress) | Max. 1000 mm         |
| Trendelenburg                     | 25°                  |
| Reverse Trendelenburg             | 25°                  |
| Lateral Tilt                      | +20°                 |
| Back Section                      | +70°                 |
| Head Section                      | +60°                 |
| Leg Section                       | 0°-90°               |
| Floor Locking System              | Center Floor Locking |

#### Salient Features

- Center Pillar • Stainless Fitting - Detachable, head, leg & Pelvic section • In built kidney bridge • FiveSectional Radio -Translucent Top.
- Base & Cover made of attractive, easy to clean Stainless Steel.
- Sophisticated mechanics provide smooth, step less articulation of the table Top for precise patient positioning. • Standby Battery Back-up with recharging circuit.





**AMRAMED**  
TECH Pvt. Ltd.

**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C



**Certified Register**

## ELECTRIC OT TABLE

**STEEL QUALITY 304 MEDICAL GRADE STAINLESS STEEL**



### Special Features

- We are providing SMPS based control box which can handle better Voltage Fluctuations, When SMPS is idle there is zero current consumption and saves power. These are as per new European and American standards where the idle power consumption should be zero and where transformer based system are outdated.
- We provide remote which is ergonomically designed to have an easy and better grip with spiral cable and it has LED backlit switch with symbolic position figures making it convenient to use even in the dark.
- We provide zero position facility i.e by pressing the single button OT Table will come to normal position.
- We provide patient reverse orientation mode when head and leg section are interchanged.
- We provide 2 memory function in control panel.
- We provide indication of battery back in the second handset.
- The Handset is designed with Back LED Indicator for better visibility.
- Head and foot section is manually operated by the means of Ratchet System.
- We provide Stainless steel Covered Base and Column Covers for easy cleaning and hygiene.
- We provide Anti Bacterial and Anti Fungal Mattress.
- The motion of table is easy with double swivel castors 100mm diameter, in which 2 castors having brak

### AMRA MATRIX ADVANCE

### Technical Data

|                                   |                      |
|-----------------------------------|----------------------|
| Length of Table                   | 1900 mm              |
| Width of Table Top                | 520 mm               |
| Minimum Height (without Mattress) | Min. 750 mm          |
| Maximum Height (without Mattress) | Max. 1000 mm         |
| Trendelenburg                     | 25°                  |
| Reverse Trendelenburg             | 25°                  |
| Lateral Tilt                      | +20°                 |
| Back Section                      | +70°                 |
| Head Section                      | +60°                 |
| Leg Section                       | 0°-90°               |
| Floor Locking System              | Center Floor Locking |

### Salient Features

- Sliding Top • Split Leg Section • Center Pillar • Stainless Fitting - Detachable, head, leg & Pelvic section • In built kidney bridge • Five Sectional Radio - Translucent Top.
- Base & Cover made of attractive, easy to clean Stainless Steel.
- Sophisticated mechanics provide smooth, step less articulation of the table Top for precise patient positioning. • Standby Battery Back-up with recharging circui



**AMRAMED**  
TECH Pvt. Ltd.

**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C



## WORKING POSITION



Height Adjustment - Down



Height Adjustment - Up



Lateral Tilt - Left



Lateral Tilt - Right



Trendelenburg



Reverse Trendelenburg



Flex Position



Gynec Position



**AMRAMED**  
TECH Pvt. Ltd.

**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C



**Certified Register**

## C-ARM COMPATIBLE ELECTRIC OT TABLE



### AMRA TMI-1201 STEEL QUALITY 304 MEDICAL GRADE STAINLESS STEEL

#### Special Features

- Eccentric Pillars
- Stainless Fitting
- Detachable, head, leg & Pelvic section
- In built kidney bridge
- Five Sectional Radio - Translucent Top.
- Base & Cover made of attractive, easy to clean Stainless Steel.
- Sophisticated mechanics provide smooth, step less articulation of the table Top for precise patient positioning.
- Standby Battery Back-up with recharging ci

#### TABLE POSITIONS

- Height Adjustment - Up & Down
- Lateral Tilt - Left & Right
- Trendelenburg & Reverse Trendelenburg
- Flex Position & Reflex Position

#### Technical Data

|                                   |              |
|-----------------------------------|--------------|
| Length of Table                   | 1900 mm      |
| Width of Table Top                | 520 mm       |
| Minimum Height (without Mattress) | Min. 750 mm  |
| Maximum Height (without Mattress) | Max. 1000 mm |
| Trendelenburg                     | 25°          |
| Reverse Trendelenburg             | 25°          |
| Lateral Tilt                      | +20°         |
| Back Section                      | +60°         |
| Head Section                      | +60°         |
| Leg Section                       | 0°-90°       |
| Floor Locking System              | Manual       |

#### OPTIONAL



Low Height



Floor Locking by Remote



Manual Override





**AMRAMED**  
TECH Pvt. Ltd.

**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C



**Certified Register**

## C-ARM COMPATIBLE ELECTRIC OT TABLE



**STEEL QUALITY 304 MEDICAL GRADE STAINLESS STEEL**



**ADDITIONAL CONTROL PANE**



**MEMORY FUNCTION**



**ZERO POSITION**

### AMRA TMI-1201 ADVANCE

#### Special Features

- We are providing first SMPS based control Box.
- Can handle better Voltage Fluctuations.
- When idle has zero current consumption and saves power, these are as per new European and American Standards where the idle power consumption should be zero and where Transformer based systems are outdated
- We are providing a Handset and Additional control panel
- Indication of Battery Back provided in the second handset.
- The Control panel has 2 memory functions.
- The second handset has also safety key to lock the complete handset
- The Handset is designed with Back LED Indicator for better visibility.
- This is the only handset with Backlight features

#### Technical Data

|                                   |              |
|-----------------------------------|--------------|
| Length of Table                   | 1900 mm      |
| Width of Table Top                | 520 mm       |
| Minimum Height (without Mattress) | Min. 750 mm  |
| Maximum Height (without Mattress) | Max. 1000 mm |
| Trendelenburg                     | 25°          |
| Reverse Trendelenburg             | 25°          |
| Lateral Tilt                      | +20°         |
| Back Section                      | +60°         |
| Head Section                      | +60°         |
| Leg Section                       | 0°-90°       |

#### TABLE POSITIONS

- Height Adjustment - Up & Down
- Lateral Tilt - Left & Right
- Trendelenburg & Reverse Trendelenburg
- Flex Position & Reflex Position



Low Height



Floor Locking by Remote



Manual Override

#### Salient Features

- Eccentric Pillars • Stainless Fitting • Detachable, head, leg & Pelvic section • In built kidney bridge • Five Sectional Radio - Translucent Top. • Base & Cover made of attractive, easy to clean Stainless Steel. • Sophisticated mechanics provide smooth, step less articulation of the table Top for precise patient positioning. • Standby Battery Back-up with recharging circuit

☎ 8275140877 📞 9822688071 ✉ info@amramedtech.com 🌐 amramedtech.com

Pune Maharashtra India

## C-ARM COMPATIBLE HYDRAULIC OT TABLE



**STEEL QUALITY 304 MEDICAL GRADE STAINLESS STEE**

### AMRA TMI-1202

#### Special Features

- Encentric Pillars
- Stainless Fitting
- Detachable, head, leg & Pelvic section
- In built kidney bridge
- Five Sectional Radio - Translucent Top.
- Base & Cover made of attractive, easy to clean Stainless Steel.
- Sophisticated mechanics provide smooth, step less articulation of the table Top for precise patient positioning.

#### Technical Data

|                                   |              |
|-----------------------------------|--------------|
| Length of Table                   | 1900 mm      |
| Width of Table Top                | 520 mm       |
| Minimum Height (without Mattress) | Min. 750 mm  |
| Maximum Height (without Mattress) | Max. 1000 mm |
| Trendelenburg                     | 25°          |
| Reverse Trendelenburg             | 25°          |
| Lateral Tilt                      | +20°         |
| Back Section                      | +70°         |
| Head Section                      | +60°         |
| Leg Section                       | 0°-90°       |
| Floor Locking System              | Manual       |

## C-ARM COMPATIBLE SEMI-ELECTRIC OT TABLE



OPTIONAL: BATTERY BACKU



**STEEL QUALITY 304 MEDICAL GRADE STAINLESS STEEL**

### AMRA TMI-1202 EL

#### Special Features

- Electric up & down by cable remote
- Eccentric Pillars
- Stainless Fitting
- Detachable, head, leg & Pelvic section
- In built kidney bridge
- FiveSectional Radio - Translucent Top.
- Base & Cover made of attractive, easy to clean Stainless Steel.
- Sophisticated mechanics provide smooth, step less articulation of the table Top for precise patient positioning.

#### Technical Data

|                                   |              |
|-----------------------------------|--------------|
| Length of Table                   | 1900 mm      |
| Width of Table Top                | 520 mm       |
| Minimum Height (without Mattress) | Min. 750 mm  |
| Maximum Height (without Mattress) | Max. 1000 mm |
| Trendelenburg                     | 25°          |
| Reverse Trendelenburg             | 25°          |
| Lateral Tilt                      | +20°         |
| Back Section                      | +70°         |
| Head Section                      | +60°         |
| Leg Section                       | 0°-90°       |
| Floor Locking System              | Manual       |





**AMRAMED**<sup>®</sup>  
TECH Pvt. Ltd.

**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C



**Certified Register**

## GENERAL SURGERY OT TABLE HYDRAULIC

### AMRA TMI-1203



HYDRAULIC UP & DOWN

## GENERAL SURGERY OT TABLE SEMI-ELECTRIC

### AMRA TMI-1203 EL



OPTIONAL: BATTERY BACKUP  
ELECTRIC UP & DOWN





## POSITION OF DILEVERY OT TABLE



### Special Features

- Center Pillar
- Stainless Fitting
- Detachable, head, leg & Pelvic section
- In built kidney bridge
- FiveSectional Radio - Translucent.
- Base & Cover made of attractive, east to clean Stainless Steel.
- Sophisticated mechanics provide smooth, step less articulation of the table Top for precise patient positioning.

### Technical Data

|                                   |              |
|-----------------------------------|--------------|
| Length of Table                   | 1900 mm      |
| Width of Table Top                | 520 mm       |
| Minimum Height (without Mattress) | Min. 750 mm  |
| Maximum Height (without Mattress) | Max. 1000 mm |
| Trendelenburg                     | 25°          |
| Reverse Trendelenburg             | 25°          |
| Lateral Tilt                      | +20°         |
| Back Section                      | +70°         |
| Head Section                      | +60°         |
| Leg Section                       | 0°-90°       |
| Floor Locking System              | Manual       |

## SLIDING TOP ELECTRIC OT TABLE



### SLIDING TOP

### STEEL QUALITY 304 MEDICAL GRADE STAINLESS STEEL

### AMRA TMI-1204

#### Special Features

This Electromatic sliding top O.T Table is suitable for all kind of surgeries, i.e. general surgery, vascular surgery, cardio, neurology, urology, gynecology, proctology, laparoscopy, traumatic surgery, plastic surgery and others. The base & column of the table is made of high quality medical stainless steel (Grade - 304). Non-reflecting surface, antibacterial & easy to clean. We provide with a top sliding of 5 inch back/ forward. Height adjustment, Lateral Tilts, Trendelenburg, and Flex - Reflex, Chair Position and Longitudinal Sliding top are realized by Remote control

#### Technical Data

|                                   |                         |
|-----------------------------------|-------------------------|
| Length of Table                   | 1900 mm                 |
| Width of Table Top                | 520 mm                  |
| Minimum Height (without Mattress) | Min. 750 mm             |
| Maximum Height (without Mattress) | Max. 1000 mm            |
| Trendelenburg                     | 25°                     |
| Reverse Trendelenburg             | 25°                     |
| Lateral Tilt                      | +20°                    |
| Back Section                      | +20°                    |
| Head Section                      | +60°                    |
| Leg Section                       | 0°-90°                  |
| Top Sliding                       | 5" Reverse & 5" Forward |

### TABLE POSITIONS

- Height Adjustment - Up & Down
- Lateral Tilt - Left & Right
- Trendelenburg & Reverse Trendelenburg
- Flex Position & Reflex Position

### Salient Features

- Eccentric Pillars • Stainless Fitting • Detachable, head, leg & Pelvic section • Four Sectional Radio - Translucent Top.
- Base & Cover made of attractive, easy to clean Stainless Steel. • Sophisticated mechanics provide smooth, step less articulation of the table Top for precise patient positioning. • Standby Battery Back-up with recharging circuit.

#### OPTIONAL



Low Height



Floor Locking by Remote



## ELECTRIC OT TABLE WITH SLIDING TOP



ADDITIONAL CONTROL PAN



MEMORY FUNCTION



ZERO POSITION

**STEEL QUALITY 304 MEDICAL GRADE STAINLESS STEEL**

### AMRA TMI-1204 ADVANCE

#### Special Features

This Electromatic sliding top O.T Table is suitable for all kind of surgeries, i.e. general surgery, vascular surgery, cardio, neurology, urology, gynecology, proctology, laparoscopy, traumatic surgery, plastic surgery and others. The base & column of the table is made of high quality medical stainless steel (Grade - 304). Non-reflecting surface, antibacterial & easy to clean. We provide with a top sliding of 5inch back/ forward. Height adjustment, Lateral Tilts, Trendelenburg, and Flex - Reflex, Chair Position and Longitudinal Sliding top are realized by Remote control

#### Technical Data

|                                   |                         |
|-----------------------------------|-------------------------|
| Length of Table                   | 1900 mm                 |
| Width of Table Top                | 520 mm                  |
| Minimum Height (without Mattress) | Min. 750 mm             |
| Maximum Height (without Mattress) | Max. 1000 mm            |
| Trendelenburg                     | 25°                     |
| Reverse Trendelenburg             | 25°                     |
| Lateral Tilt                      | +20°                    |
| Back Section                      | +20°                    |
| Head Section                      | +60°                    |
| Leg Section                       | 0°-90°                  |
| Top Sliding                       | 5" Reverse & 5" Forward |

#### TABLE POSITIONS

- Height Adjustment - Up & Down
- Lateral Tilt - Left & Right
- Trendelenburg & Reverse Trendelenburg
- Flex Position & Reflex Position

#### Salient Features

- Eccentric Pillars • Stainless Fitting • Detachable, head, leg & Pelvic section • Four Sectional Radio - Translucent Top.
- Base & Cover made of attractive, easy to clean Stainless Steel. • Sophisticated mechanics provide smooth, step less articulation of the table Top for precise patient positioning. • Standby Battery Back-up with recharging circuit.

#### OPTIONAL



Low Height



Floor Locking by Remote

## NEURO ATTACHMENT



**STEEL QUALITY 304 MEDICAL GRADE STAINLESS STEEL**

### AMRA TMI-1205

We are ranked as one of the leading manufacturers, suppliers, exporters, traders and distributors of excellent quality Neuro Surgery Attachments. These neuro surgery attachments provide neuro surgeons stable and excellent support as they perform different medical procedures. These neuro surgery attachments are very durable made from excellent quality stainless steel and other materials. We provide ophthalmologic procedures using the narrow headrest, neurology procedures in prone position and supine position (Sugita six pin fixator and Horse Shoe).



Horse Shoe Position



Sugita Head Fixator for  
Supine position & Skull Pins



**AMRAMED**  
TECH Pvt. Ltd.

**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C



**Certified Register**

## ALBEE'S TYPE ORTHO ATTACHMENT



### AMRA TMI-1206 ALBEE'S TYPE

#### Special Features

Ortho Attachment allows the table to be used for the lower extremity procedures. Available at competitive prices our ortho attachment is purely made of stainless steel which resist corrosion.

#### Salient Features

- The Orthopedic Attachment can be used with any operating table.
- In addition, a pelvic support plate is supported on top of the frame and is formed with an asymmetrical construction.
- Stirring Instruments allows flexion/raising of legs & stirring/string ball joints locked from foot end.
- Two traction devices on sliding spot with ball joint allowing flexions & rotations of foot with traction shoe with metallic sole.
- Padded orthopedic knee crutches/provides support to the legs of the patient.



Patient positioning for Neck of femur nailing traction (D.H.S)



Femur treatment in lateral position



Tibia extension position





**AMRAMED**  
TECH Pvt. Ltd.

**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C



**Certified Register**

## HANGING TYPE ORTHOPEDIC ATTACHMENT



### AMRA TMI-1206 HANGING TYPE

#### Special Features

Ortho Attachment allows the table to be used for the lower extremity procedures. Available at competitive prices our ortho attachment is purely made of stainless steel which resist corrosion.

#### Salient Features

- The Orthopedic Attachment can be used with any operating table.
- In addition, a pelvic support plate is supported on top of the frame and is formed with an asymmetrical construction.
- Stirring Instruments allows flexion/raising of legs & stirring/string ball joints locked from foot end.
- Two traction devices on sliding spot with ball joint allowing flexions & rotations of foot with traction shoe with metallic sole.
- Padded orthopedic knee crutches/provides support to the legs of the patient.



**AMRAMED**  
TECH Pvt. Ltd.

**Wholesale Supplier**


A to Z Hospital Equipment and  
Surgical Products B to C



**Certified Register**

## OPHTHALMIC ELECTRIC OT TABLE

**STEEL QUALITY 304 MEDICAL GRADE STAINLESS STEEL**

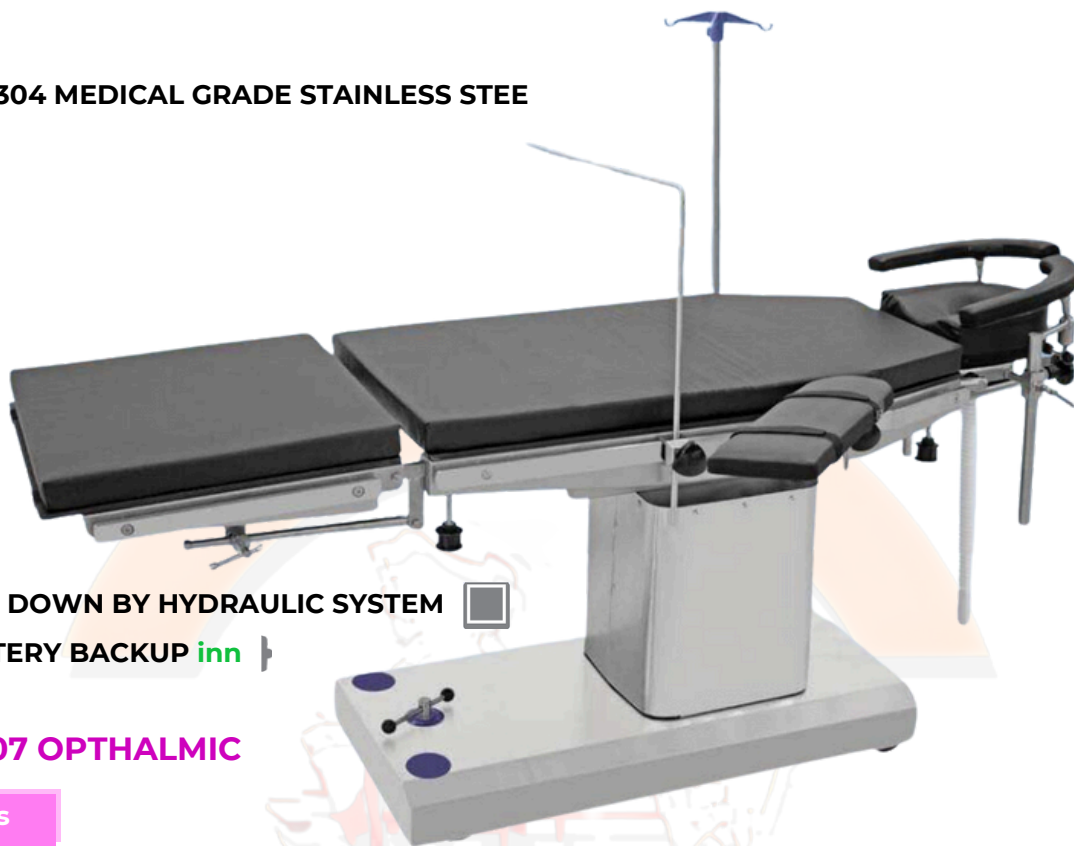

**OPTIONAL: UP & DOWN BY HYDRAULIC SYSTEM**



**OPTIONAL: BATTERY BACKUP**



### AMRA TMI-1207 OPTHALMIC

#### Special Features

This Ophthalmic Table is specially designed for ophthalmic surgery. Three sectional table is provided with drainage system in the head section. The table has sturdy base, column made up of stainless steel, antibacterial and is easy to clean. Electromotive height with the help of lightweight handheld control

#### Standard Accessories

- Anesthetic Screen
- IV Pole
- Hand Rest
- Mattres

#### Technical Data

|                                   |              |
|-----------------------------------|--------------|
| Length of Table                   | 1900 mm      |
| Width of Table Top                | 520 mm       |
| Minimum Height (without Mattress) | Min. 750 mm  |
| Maximum Height (without Mattress) | Max. 1000 mm |
| Leg Section                       | 0°-90°       |

#### Salient Features

• Stainless Fitting • Detachable, head & leg section • Three Sectional Radio Translucent Top. • Base & Cover made of attractive, easy to clean Stainless Steel. • Sophisticated mechanics provide smooth, step less articulation of the table Top for precise patient positioning



**Foot Switch (Optional)**



**Instrument Tray (Optional)**



**AMRAMED**  
TECH Pvt. Ltd.

**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C



**Certified Register**

## ELECTRIC GYNECOLOGICAL OBSTETRIC OT TABLE



Examination &  
Birthing Positio

### Colour options



**OPTIONAL: BATTERY BACKUP**

**STEEL QUALITY 304 MEDICAL GRADE STAINLESS STEEL**

### AMRA TMI-1208

#### Special Features

Obstetrics operating room table perfect for all kinds of gynecological operations. The cover for base and column is made of 304 stainless steel, good performance on anti-rust and easy to clean.

#### Salient Features

- Functional and versatile
- Flexibility and precision
- Hygienic design
- Easy cleaning • User-friendly

#### Standard Accessories

- Arm Rest - 2 pcs.
- Birth Handle - 2 pcs.
- Mattress
- GynaeTray



Gynae Tray



Birthing Handle

#### Technical Data

|                                   |             |
|-----------------------------------|-------------|
| Length of Table                   | 1900 mm     |
| Width of Table Top                | 520 mm      |
| Minimum Height (without Mattress) | Min. 750 mm |
| Maximum Height (without Mattress) | Max. 950 mm |
| Trendelenburg                     | ±20°        |
| Back Section                      | ±70°        |
| Floor Locking System              | Manual      |



Foot Switch (Optional)



**AMRAMED**  
TECH Pvt. Ltd.

**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C



**Certified Register**

## DERMA CHAIR

OPTIONAL: BATTERY BACKUP

**STEEL QUALITY 304 MEDICAL GRADE STAINLESS STEEL**

**AMRA TMI-1209**

### Special Features

- Luxurious Styling
- Multi-position section couch
- Suitable for all Dermatology, Cosmetology & Laser Surgery.
- Stable and professional design, unlike Saloon Chairs.
- Three independent motors - for Height, Leg Section & Back Section
- Elegant upholstery and all english colour

**Colour options**



**Foot Switch (Optional)**

### Technical Data

|           |               |
|-----------|---------------|
| Power     | 230 V 50      |
| Size (mm) | 2000x610x700  |
| Weight    | 80 Kgs        |
| Motors    | 3 independent |
| Remote    | Provided      |

**Height Adjustment - Up / Down**



## ELECTRIC ANIMAL OPERATING TABLE



**STEEL QUALITY 304 MEDICAL GRADE STAINLESS STEEL**

**OPTIONAL: UP & DOWN BY HYDRAULIC SYSTEM**



**OPTIONAL: BATTERY BACKUP** mm

**AMRA TMI-1301**

### Special Features

- Electric up and down function.
- It is suitable for small animal's operation, such as sheep, dog etc.
- It is reasonably structured, tight, high performance, easy operation and has very strong practicability.
- Fully made of Premium Grade 304 Stainless Steel, it is heat-resistant, corrosion protective and rustproof.
- The surface of the table-top is V shaped and is smooth and easy to clean and disinfect, thus, it can avoid cross infection
- Equipped with fluid collecting through system.




**Foot Switch (Optional)**

### Technical Data

|                 |                   |
|-----------------|-------------------|
| Length of Table | 1220mm (48")      |
| Width of Table  | 610mm (24")       |
| Minimum Height  | Min. 660 mm (26") |
| Maximum Height  | Max. 915 mm (36") |

## STANDARD ACCESSORIES

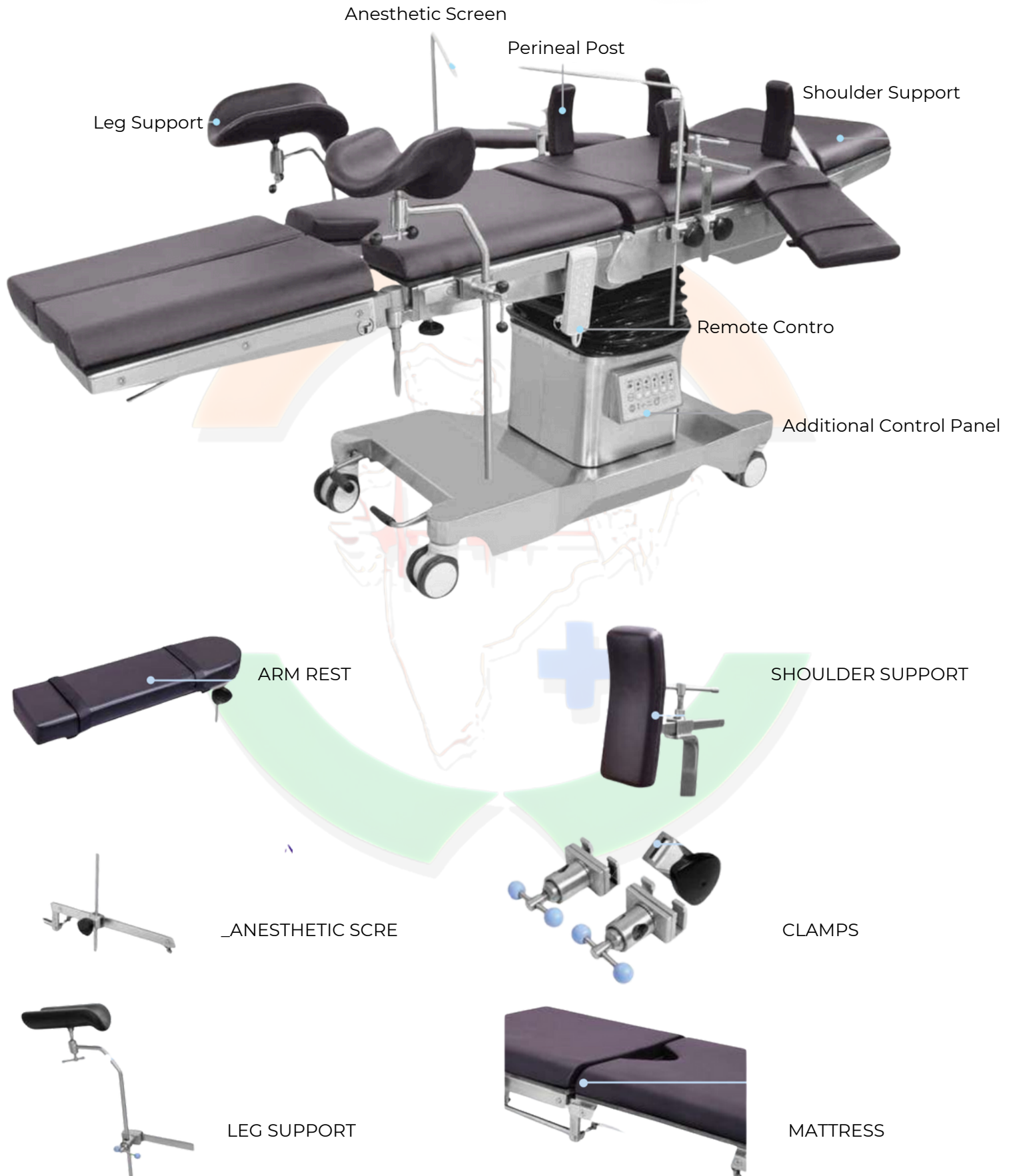



**AMRAMED**  
TECH Pvt. Ltd.

**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C



**Certified Register**

## OPTIONAL ATTACHMENT ACCESSORIES



**SPLIT LEG SECTION**  
OA-1



**RESTRAINT BODY STRAP**  
OA-2



**BOLSTER PILLOW**  
OA-3



**INSTRUMENT TRAY**  
OA-4




**FOOT SUPPORT**  
OA-5



**CERVICAL IN SUPINE /  
LAMINECTOMY FRAME**  
OA-6



**HUMERUS SUPPORT**  
OA-7

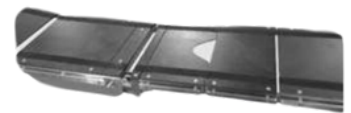

**TABLE LENGTH EXTENSION**  
OA-8



**STEINMANN PIN  
CLAMP HOLDER**  
OA-9



**SPINE SURGERY USING  
CUSHION SPINE FRAME**  
OA-10



**ARM TRACTION DEVICE  
/ PERINEAL POS**  
OA-11



**UROLOGY I GYNAE TRAY**  
OA-12



**HORSE SHOE**  
OA-13



**KNEE SUPPORT**  
OA-14




**FIXATOR FOR  
SITTING POSITION**  
OA-15



**SHOULDER SURGE**  
OA-16



**HAND SURGERY  
TABLE**  
OA-17

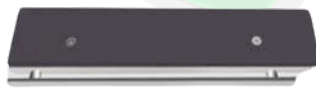

**WIDTH EXTENSION**  
OA-18



**ARTHROSCOPY  
ATTACHMENT**  
OA-19



**WRISTLET**  
OA-20



**X-RAY TOP**  
OA-21



**ORTHOPAEDIC  
ATTACHMENT STAN**  
OA-22



**ACCESSORIES STAN**  
OA-23