


Linked Page



Certified Register

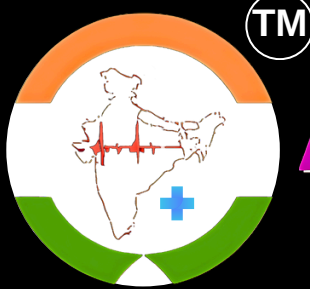

**AMRAMED**  
**TECH Pvt. Ltd.**

**Wholesale Supplier**  
**Operation Table Economical**  
**A to Z Hospital Equipment and**  
**Surgical Products B to C**



 **8275140877**

 **9822688071**

 **info@amramedtech.com**



**Sr No.29/1/2, Adarsh Nagar,  
Uruli Devachi, Hadapsar, Pune**  
**Branches :- All Maharashtra**



**amramedtech.com**

## SURGICAL OT TABLES

It is single action O.T. table with improved hydraulic system provides smooth uniform movement and the speed of lowering the table is permanently set. All the procedures positions can be attained through a crank. The detachable leg section provides easy lithotomy positioning. Head and leg section are interchangeable. Heavy mild base and column are covered with stainless steel.



### AMRA-3000

HYDRAULIC	202 S S
HYDRAULIC	304 S S
SEMI-ELECTRIC	202 S S
SEMI-ELECTRIC	304 S S

#### Technical Parameters

Table Length	<b>1900 mm</b>
Width	<b>500 mm</b>
Minimum Height	<b>750 mm</b>
Maximum Height	<b>1000 mm</b>
Kidney Elevator	<b>120 mm</b>
Lateral Tilt	<b>+20°/ —20°</b>
Trendelenburg	<b>+25°/ —25°</b>
Back Section	<b>+90°/ —40°</b>
Head & Leg Section	<b>+90°/ —80°</b>
Load Capacity	<b>180 KG</b>

### Features :

- Hydraulic and Mechanical Operated.
- Centric pillars, Stainless fitting.
- Leg & Head section are interchangeable.
- Detachable head & leg section.
- Inbuilt kidney bridge.
- Base and cover made of attractive, easy to clean stainless steel.



## (C- Arm Compatible)

C- Arm Compatible operating table is fully manually operated for all kind of surgical operations and treatments ( Orthopedics, Cardio, Neuro and Endoscopy), characterized with ergonomy and esthetic design. It is hydraulic operated in elevation and gear-mechanism operated in all types of different positions during surgical procedures.



### AMRA-5000

HYDRAULIC	202 S S
HYDRAULIC	304 S S
SEMI-ELECTRIC	202 S S
SEMI-ELECTRIC	304 S S

**\*ELECTRIC UP & DOWN FEATURE IS ALSO AVAILABLE IN THIS MODEL**

#### Technical Parameters

Table Length	<b>1950 mm</b>
Width	<b>525 mm</b>
Minimum Height	<b>750 mm</b>
Maximum Height	<b>1050 mm</b>
Kidney Elevator	<b>120 mm</b>
Lateral Tilt	<b>+20° / - 20°</b>
Trendelenburg	<b>+25° / - 25°</b>
Back Section	<b>+90° / - 40°</b>
Head & Leg Section	<b>+90° / - 80°</b>
Load Capacity	<b>200 KG</b>

### Features :


- Fully electrically operated table.
- Six sections, radio-translucent table with eccentric pillar.
- Interchangeable head and leg sections (manually)
- All positions operated through cable remote.
- Inbuilt kidney bridge.
- Support on a heavy base column covered with stainless steel sheet.
- Full body imaging without special accessories.



**Knee Crutches (1 pair)**




**Padded Shoulder Supports (1 Pair)**

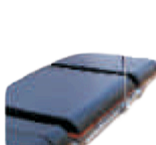


**Anaesthetic Frame (1Pcs)**



**Intravenous Armboard (1 pair)**

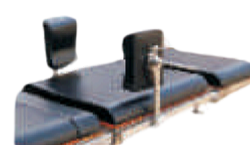

**Padded Side Supports (1Pair)**

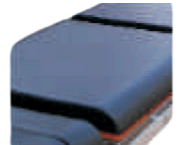

**One Mattress Set Thickness 50mm**



## ELECTRO TABLE

An Innovating arrange of C- Arm compatible multi utility electric operating table equipped with intuitive intelligence & Designed with minute articulations for ergonomic convenience imaging capability & attachable accessories, to enhance comfort & precision for the surgical team



### AMRA-8000

FULLY-ELECTRIC 202 S S

FULLY-ELECTRIC 304 S S

#### Technical Parameters

Table Length	1950 mm
Width	525 mm
Minimum Height	750 mm
Maximum Height	1050 mm
Kidney Elevator	120 mm
Lateral Tilt	+20° / - 20°
Trendelenburg	+25° / - 25°
Back Section	+90° / - 40°
Head & Leg Section	+90° / - 80°
Load Capacity	200kg
Battery Backup	2-3 Hrs. (Optional)
Mains Voltage	AC 220V + 10% 50Hz

### Features :

- Fully electrically operated table.
- Six sections, radio-translucent table with eccentric pillar.
- Interchangeable head and leg sections (manually)
- All positions operated through cable remote.
- Inbuilt kidney bridge.
- Support on a heavy base column covered with stainless steel sheet.
- Full body imaging without special accessories.



## ADVANCE ELECTRO TABLE

**Perfect functionality and versatility.** The ADVANCE ELECTRO TABLE offers effective solutions for the widest variety of requirements in the OT. Patients benefit from superb positioning comfort and easily accessed with corded handset. **Universal Design.** Developed with numerous new safety features, the new table is straightforward to use in all surgical disciplines. A wide range of patient positions makes this operating table universal. **Maximum Performance for minimum investment.**



### AMRA-10000

F ADV ELECTRIC	202 S S
WITHOUT IMPORTED WHEELS	
F ADV ELECTRIC	304 S S
WITH IMPORTED WHEELS	

#### Technical Parameters

Table Length	<b>1950 mm</b>
Width	<b>525 mm</b>
Minimum Height	<b>750 mm</b>
Maximum Height	<b>1050 mm</b>
Kidney Elevator	<b>120 mm</b>
Lateral Tilt	<b>+20° / - 20°</b>
Trendelenburg	<b>+25° / - 25°</b>
Back Section	<b>+90° / - 40°</b>
Head & Leg Section	<b>+90° / - 80°</b>
Load Capacity	<b>250 KG</b>
Battery Backup	<b>2-3 Hrs. (Optional)</b>
Mains Voltage	<b>AC 220V + 10% 50Hz</b>

### Features :

- SS-304
- Surgeons Panel
- Electric (remote operated)
- C Arm Compatible
- Zero Positioning
- Casters lockable with single action break
- Handset is providing Dual Memory Function which can be set and reset.



**\* Top Slide Features also Available**

## GYNAECOLOGY & OBSTETRICS OT TABLE

In terms of quality, reliability and functionality, this Operating Table is being acknowledged as the ideal solution for Gynaecological diagnosis and examination, child-birth, cesarian operation and other gynaecological surgical operations. Due to its modular design, it can be flexibly adapted to your needs and configured to suit your exact surgical requirements.

### Modular Construction


provides flexibility and economy

### Ergonomic Functionality

for an efficient and safe application

### High Stability

High load bearing capacity and a long service life

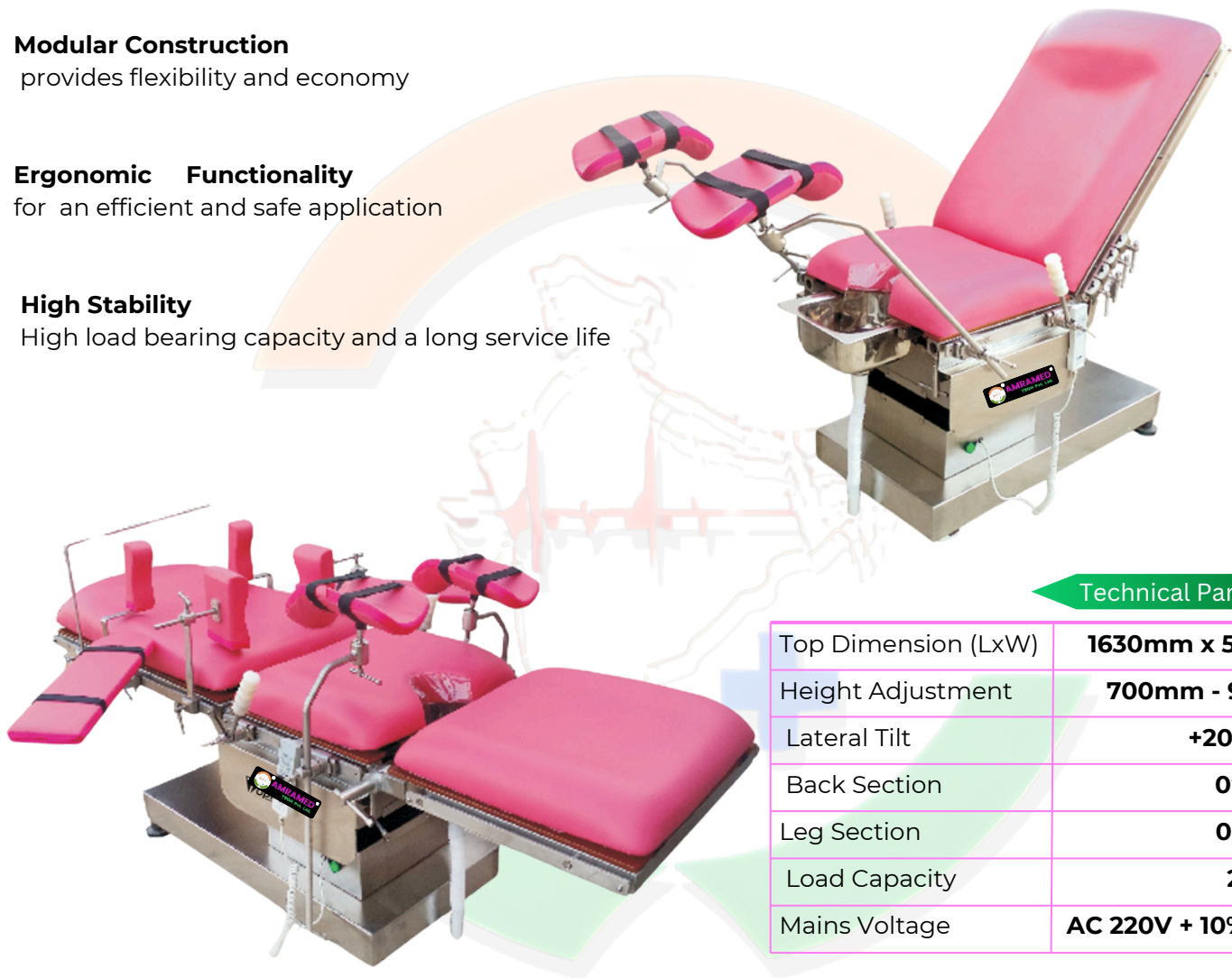

### Technical Parameters

Top Dimension (LxW)	<b>1630mm x 550 mm</b>
Height Adjustment	<b>700mm - 950mm</b>
Lateral Tilt	<b>+20°/ —20°</b>
Back Section	<b>0°/ —80°</b>
Leg Section	<b>0°/ —90°</b>
Load Capacity	<b>200 KG</b>
Mains Voltage	<b>AC 220V + 10% 50Hz</b>

### Standard Accessories :-

- Knee Crutch with Straps : One Pair
- Arm Board with Mattress : One Pair
- Padded Side & Shoulder Support : One Pair
- Hand Grip with Clamp : One Pair
- SS Douche Tray : One Pc.
- 

For special order we can change technical parameters,  
not changing the safety features of the table.

**AMRA-GO-H**  
**AMRA-GO-E**

## OPHTHALMIC TABLE

This O.T table is specially designed for Ophthalmic surgery . Three sectional table is provided with drainage system in the head section. The table has sturdy base, column made up of Stainless steel, antibacterial and is easy to clean. Electromotive height and trendelenburg positions with the help of lightweight handheld control unit. Easy detachable leg and head section.



AMRA-EYE-H  
AMRA-EYE-SE  
AMRA-EYE-FE

### Technical Parameters

Length Of Table top	1950mm
Width of table top	500mm
Minimum Height	700mm
Maximum Height	950mm
Trendelenburg	+20° /- 30°
Load Capacity	180 KG



Maximum Height  
950 mm



Minimum Height  
700 mm



Trendelenburg  
20°



Reverse Trendelenburg  
30°

**\*Foot Remote Optional.**

**“For special order we can change technical parameters,  
not changing the safety features of the table.”**



## DERMATOLOGY AND HAIR TRANSPLANT CHAIR

The Dermatology chair we supply in the market is used in various Hair Transplant, Dermatology and Cosmetology Centres. This Hair Transplant Chair/Derma Bed comes with motorized positioning control that is intended for medical purpose and that can be adjusted to various positions through a hand-operated remote. The table has an aesthetically beautiful bodyline and is fully covered with non-corrosive, mat finish SS 304 Grade Stainless Steel for easy cleaning and hygiene.

### UNIQUE FEATURES:

- Fully Automatic Remote controlled chair with noiseless imported motor system.
- The chair is capable of handling a load of upto 250kg.
- Three Actuators system allows independent Up/Down, Back Rest & Foot Rest movements.
- Elegant Upholstery with custom colours.
- PU-Mattress, Better than foam - mattress.
- Imported Castor wheels, easy to move.
- Stainless Steel material (Grade 202).
- Special Head Rest for Hair Transplant with angle adjustment.

### Technical Parameters

Loading Capacity	200kg
Motor	3 independent Motors
Colour	Blue, Ocean Green, Gray, White, Black & Pink.
Weight	125kg
Ambient Temp.	5° C to 45° C
Size(in mm)	1950 X 625 X 700
Power	230 V 50 Hz



### POSITIONS :



Maximum Height  
800 mm



Minimum Height  
650 mm



Foot Rest Movement  
75°



Back Rest  
80°



Chair Position

### AMRA DHTC 5000

AMRA-DHTC-M

AMRA-DHTC-H

AMRA-DHTC-E2

AMRA-DHTC-E3

**\*Foot Remote Optional**

**\*For special order we can change technical parameters.  
not changing the safety features of the chair.**

## DERMATOLOGY AND HAIR TRANSPLANT CHAIR


The Dermatology chair we supply in the market is used in various Hair Transplant, Dermatology and Cosmetology Centres. This Hair Transplant Chair/Derma Bed comes with motorized positioning control that is intended for medical purpose and that can be adjusted to various positions through a hand-operated remote. It is very much stable and is made of complete stainless steel.

### Unique Features :

1. Fully Automatic Remote controlled chair with noiseless imported motor system.
2. Four Actuators system allows independent Up/Down , Back Rest, Trendelenburg & Foot Rest movements.
3. Trendelenburg Position easy to use for long duration procedures with chair position.
4. PU-Mattress ,Better than form - mattress.
5. Imported Castor wheels, easy to move.
6. Stainless Steel material (Grade 304).
7. Special Head Rest for Hair Transplant with angle adjustment

### Technical Parameters

Loading Capacity	200kg
Motor	4 independent Motors
Colour	Blue, Ocean Green, Gray, White, Black & Pink.
Weight	150kg
Ambient Temp.	_5° C to 45° C
Size(in mm)	1950 X 625 X 700
Power	230 V 50 Hz

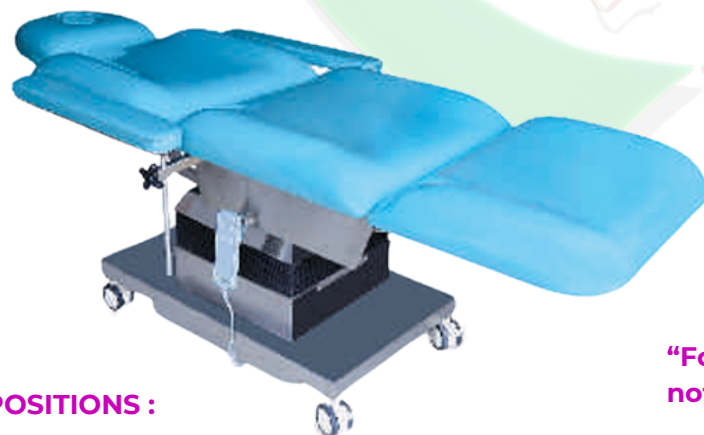

## AMRA PREMIUM DHTC 10000

AMRA-DHTC-M

AMRA-DHTC-H

AMRA-DHTC-E2

AMRA-DHTC-E3

**"For special order we can change technical parameters, not changing the safety features of the chair."**

### POSITIONS :



**Maximum Height**  
800mm



**Minimum Height**  
650mm



**Foot Rest Movement**  
75



**Back Rest**  
80




**Chair Position**



**Trendelenburg with Chair Position**

## FULLY BARIATRIC OT TABLE WITH DIVIDED LEG SECTION



It is specially designed for Bariatric Surgery. Being Robust and reliable the O.T table can bear load of heavy patient of up to 300 kgs, Its split leg section (0° to 180°) provides the surgeon wide and easy access between legs of the patient. Its steep reverse trendelenburg position surgeon can easily sit during long surgery. Supported with various accessories it offers a complete solution for range of surgeries on bariatric patients.

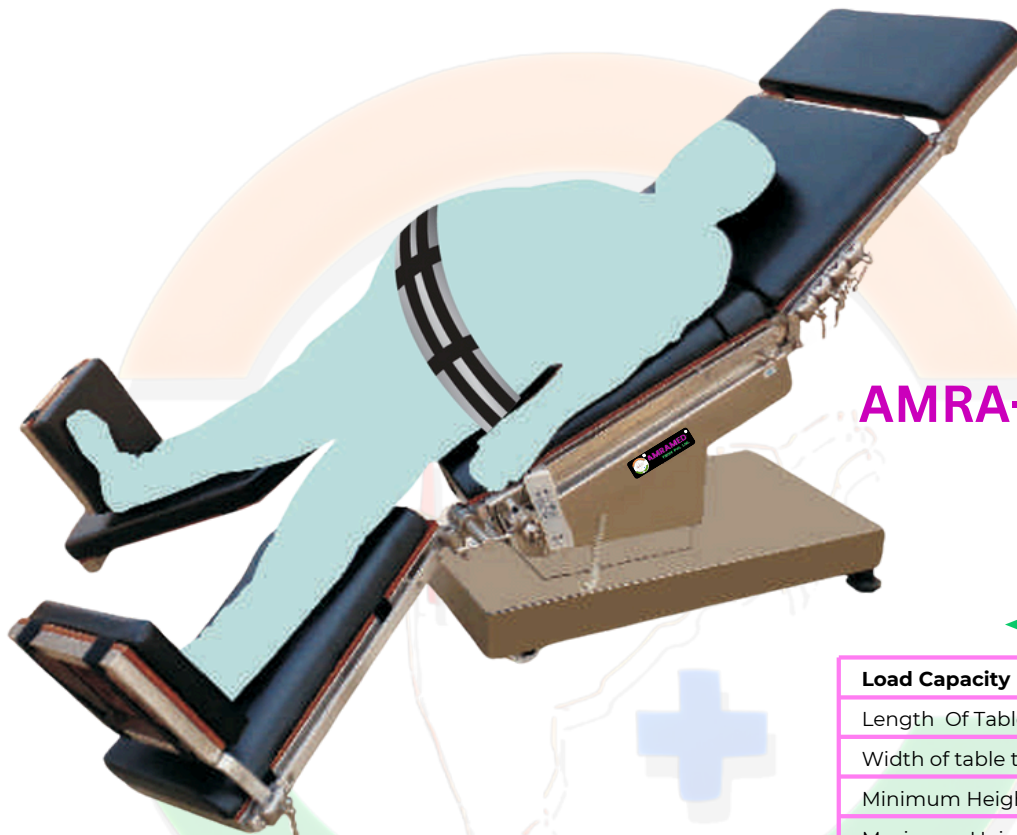

### AMRA-BARIATRIC

#### Technical Parameters

<b>Load Capacity</b>	<b>300kg</b>
Length Of Table top	2150mm
Width of table top	600mm
Minimum Height	600mm
Maximum Height	800mm
Trendelenburg	+20° / -35°
Side Tilt	20°
Leg Section (Up / Down)	+90° / -90°
Back Section	+75° / -15°

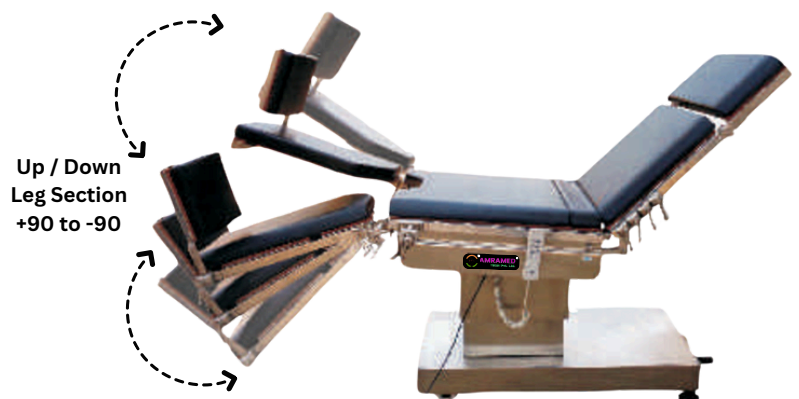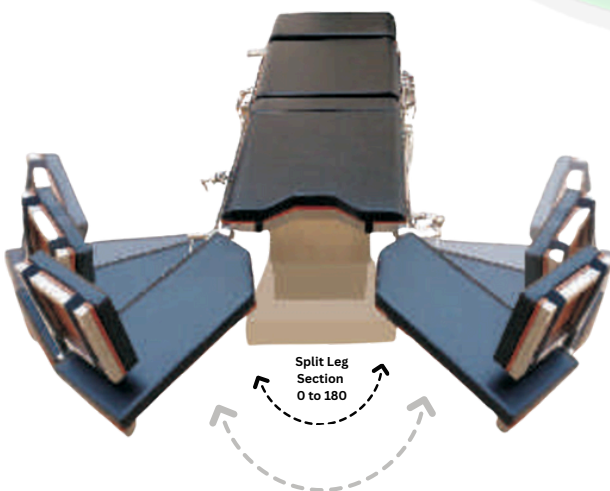





**AMRAMED**<sup>®</sup>  
TECH Pvt. Ltd.

**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C



**Certified Register**

## ELECTRIC 7 HYDRAULIC SURGEON'S CHAIR



- Surgeon while sitting can raise/lower the height.
- Minimum height 23" and maximum 32".
- Detachable Arm Rests.
- Adjustable back rest has vertical and horizontal movements.
- Full swiveling castors for easy mobility.
- Castors with brake optional.
- Complete Stainless Steel parts with upholstered seat.


## OPTIONAL ACCESSORIES

### Orthopaedic Attachment



1. Pelvic Rest
2. Pelvic Support (Radiolucent)
3. Inner thigh Support
4. Thigh Support
5. 1 Pair Leather Shoes
6. Tibia attachment
7. Hand Surgery Table (Radiolucent)
8. Hand Traction Device
9. Steinman Pin Holding Clamp
10. Special Head Rest for Shoulder Surgery

### Neuro Attachment

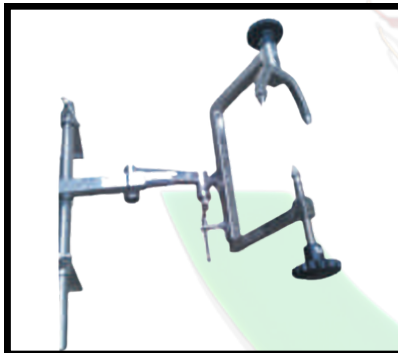

Neuro Sujita Head Fixator for  
Sitting / Prone Position



Head Rest device for prone, Spine & Lateral Position



### Accessories Trolley



### Spine Attachment



Spinal Bridge Frame



**AMRAMED**  
TECH Pvt. Ltd.


**Wholesale Supplier**

A to Z Hospital Equipment and  
Surgical Products B to C




**Certified Register**





## ELECTRIC TABLE C-ARM (Measurement)

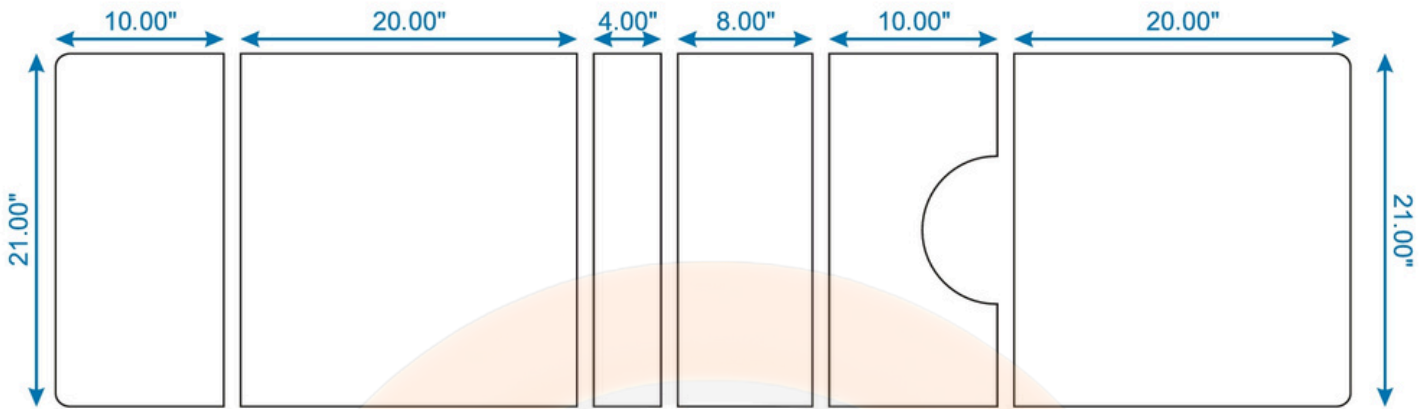

## HYDRAULIC TABLE C-ARM (Measurement)

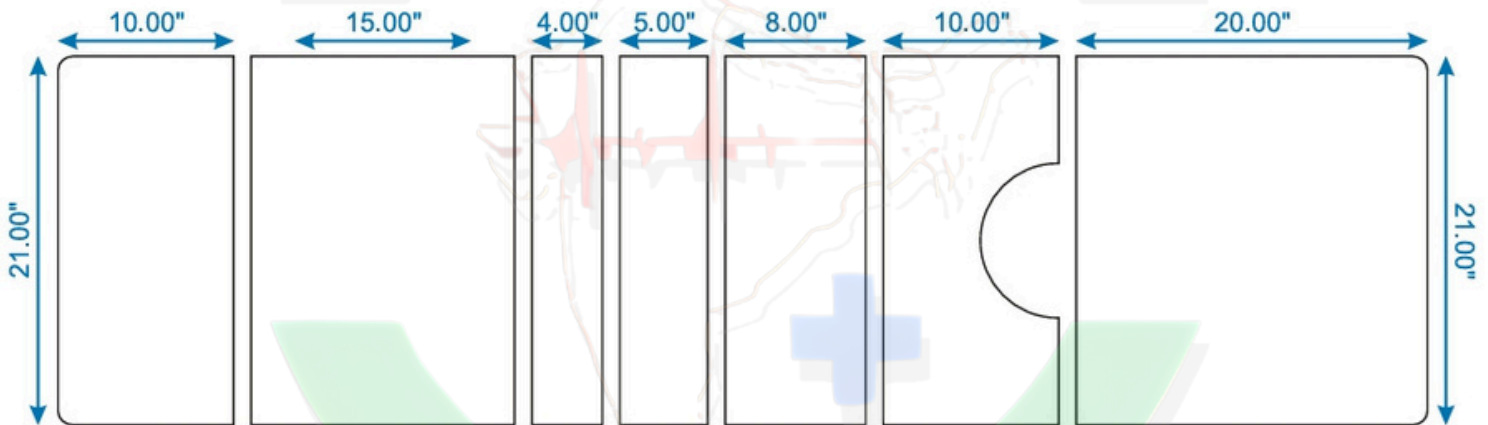

## VARIOUS POSITIONS

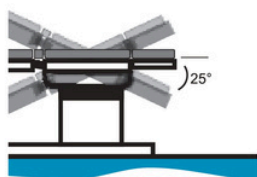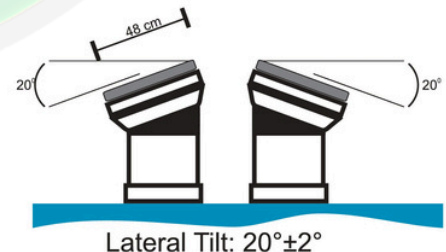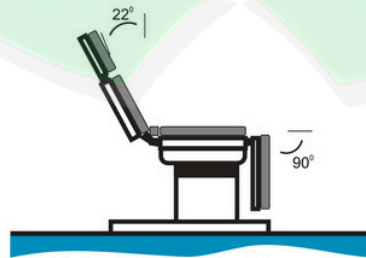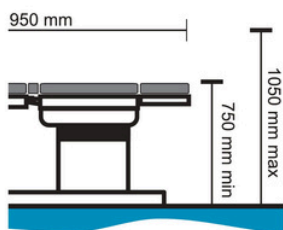

Telenburg: 25°±2°



Flex



Reflex